



WATO EX-35

Anaesthesia Machine

exceptional performance



mindray
healthcare within reach

Since 2006 Mindray has successfully installed over 20,000 anaesthetic machines to satisfied customers all over the world.

Mindray takes immense pride knowing that every few seconds a clinician, somewhere in the world, will be switching on a Mindray anaesthesia machine and using it with complete confidence and satisfaction.

Over the last decade Mindray has continued to work closely with clinicians, across the world, to recognise and understand the clinical challenges encountered everyday and overcome them with new innovative and intuitive solutions. With this in mind, Mindray is now pleased to offer the new improved WATO EX-35 with advanced features, functionality and fresh new look.



WATO EX-35 ^{NEW}

exceptional performance



New Touch-screen Control

Featuring a new 10.4" touch-screen with quick access parameter controls and flat menu design. The new WATO EX-35 UI enables the user to make adjustments to vent mode and monitor settings easier than ever before at this product level.

The new touch responsive Smart Alarm system provides real-time detailed information that quickly and clearly alerts the clinician to adverse situations highlighting the need for early interaction and resolution.

Enhanced Performance

Employing the latest technologies the new WATO EX-35 supports an extensive range of features and functionality, ensuring all anaesthetic procedures are delivered safely, efficiently and effectively.

Enhanced features include, but are not limited to:

Integrated Plug-and-play AG Module (auto ID of 5 agent and N₂O)

Integrated multi-mode ventilator offering controlled and support modes of ventilation

The new and unique Mindray Ventilator Drive Gas Auto Switch function which guarantees the patient is ventilated without interruption at all times...

Better Usability

Designed with the user in mind.

The new EX-35 has been developed to make your life easier. A feature rich, ergonomically designed anaesthesia system that can be positioned with ease and used for prolonged periods without experiencing unnecessary stress or fatigue.



New Touch-screen Control

Touch-screen, Easier Settings

The new WATO EX-35 employs a 10.4" high resolution touch-screen displaying a maximum of 3 ventilation waveforms (EtCO₂, Pressure, Flow & Volume), Spirometry Loop and full gas monitoring at the same time enabling direct access to all parameter controls and settings. The new tabular, flat level, menu structure with quick touch control is a global first for devices at this market level.

Touchable Smart Alarm Management

The new Smart Alarm management system provides real-time detailed information that quickly and clearly alerts the clinician to adverse situations and supports possible early resolution.

The comprehensive alarm log includes a user configurable alarm log showing both:

- Physical alarm information with displayed parameter snapshot and
- Technical alarms which use a diagnostic software offering possible causes for alarms that may assist in trouble-shooting in the immediate situation.



Better Usability

With conveniently located handles and smooth running castors the new WATO EX-35 can easily be positioned effortlessly within the workplace.

The large easy-clean work surface provides plenty of space for note writing, placement of drugs and intubation equipment and an ergonomically designed footrest offers the user total comfort throughout the longest of cases.





Enhanced Performance

Innovative P-n-P Gas Monitoring

Our improved Plug and Play single slot CO2 module or double slot anaesthetic agent module provides comprehensive, visual, breath-by-breath analysis for FiO2, EtO2, CO2, N2O and auto detection of 5 anaesthetic agents. What's more the modules can be moved between functionally similar Mindray monitors and anaesthetic devices further enhancing device configuration and flexibility to suit your needs.


Driven Gas Auto-switch

Safety is at the heart of all Mindray's devices and the new WATO EX-35 is no exception. The all new ventilator gas Auto-switch function ensures the patient is mechanically ventilated at all times even if the primary ventilator drive gas supply fails.



High Performance Ventilator

The new WATO EX-35 supports a full range of controlled and support modes of ventilation and with dynamic tidal volume compensation, from as low as Tve 20ml, the new WATO EX-35 is suitable for all patient categories from neonates to adult.

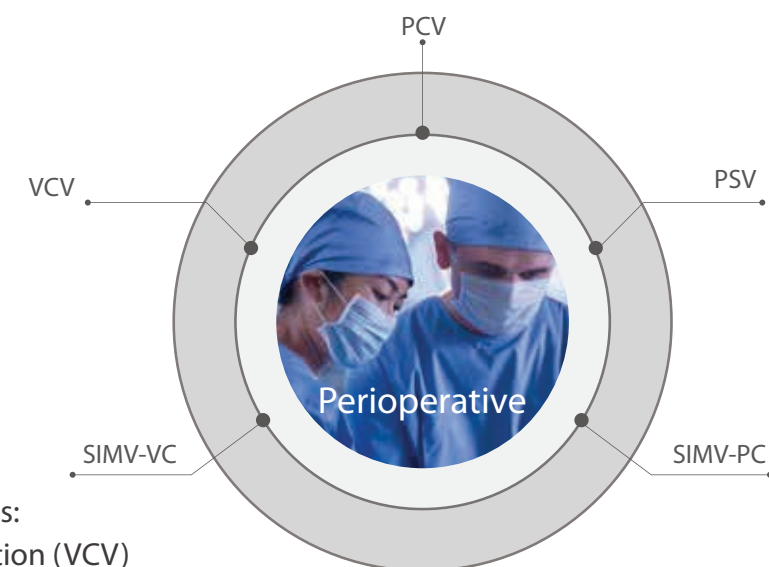

Available ventilation modes:

- Volume controlled Ventilation (VCV)
- Pressure Controlled Ventilation (PCV)
- Pressure Support Ventilation + apnea ventilation (PSV + apnea ventilation)
- Synchronised Intermittent Mandatory Ventilation – VC (SIMV -VC)
- Synchronised Intermittent Mandatory Ventilation – PC (SIMV -PC)

The materials and information provided via this brochure are distributed internationally, but not all the products or services found on this brochure are necessarily available in your country or area. Please contact local office or agents to make sure of the availability of the relevant products or services.

Mindray Medical Netherlands B.V.

Drs.W.van Royenstraat 8, 3871 AN Hoelvelaken, The Netherlands
Voice: (31-33) 254-4911
Fax: (31-33) 253-4280

Mindray Medical Italy S.r.l.

Via Leonardo DA Vinci, 158, 20090 Trezzano Sul Naviglio (MI), Italy
Voice: (39-02) 5737-401
Fax: (39-02) 5523-2018

Mindray Medical España S.L.

Avenida Manteras 38, Edificio B, Bajo B008, 28050-Madrid Spain
Voice: (34-91) 392-3754
Fax: (34-91) 392-3756

Mindray Medical Germany GmbH

Goebelstrasse 21, 3. OG, 64293 Darmstadt
Voice: (49 6151) 3910 0
Fax: (49 6151) 3910 300

Mindray (UK) Limited

Mindray House, Kingfisher Way, Hinchingsbrooke Business Park,
Huntingdon, Cambs, PE29 6FN
Voice: (44) 1480-416-840
Fax: (44) 1480-436-588

Mindray Medical France SARL

Europarc Creteil 1 allée des cerisiers 94000 CRETEIL, France
Voice: (33-1) 4513-9150
Fax: (33-1) 4513-9151

Mindray Medical Sweden AB

Rissneleden 136, 174 57 Sundbyberg, Sweden
Voice: (46) 8 555 54 121

Mindray Medical Rus Co.Ltd

2 Zvenigorodskaya Street, 13, bld 41, 123022 Moscow, Russia
Voice: (7-499) 553-6036
Fax: (7-499) 553-6039

Mindray Medical Colombia S.A.S

Room No.702, Torre UNIKA,Carrera 9 No.77-67, Bogota, Colombia
Voice: (57-1) 3130892 / 3210916

Mindray do Brasil – Comércio e Distribuição de Equipamentos Médicos Ltda.

Rua Tavares Bastos, 329, Perdizes, São Paulo - SP,
CEP 05012-020
Voice: (55-11) 3124-8026
Fax: (55-11) 3078-8035

Mindray Medical Mexico S. de R.L. de C.V.

Félix Parra # 175 Colonia. San José Insurgentes Sur Delegación Benito Juárez 03900 México, D. F.
Voice: (52-55) 5661-9450 / (52-55) 5662-6620
Fax: (52-55) 5662-6597

PT.Mindray Medical Indonesia

The East 11th Floor, Unit 08 Jl. Dr. Ide Anak Agung Gde Agung Kav. E 3.2 No.1 Jakarta 12950
Voice: (6221) 2902-7280
Fax: (6221) 2902-7283

Mindray Medical Egypt

Ground Floor -11 El Imam Metwally El Sharawy Street,
Sheraton Heliopolis - El Nozha, Cairo, Egypt
Voice: (2) 02 22691047
Fax: (2) 02 22691046

Mindray Medical Australia

ADD: 19 Harker St, Burwood VIC 3125, Australia
Voice: (61-3) 9830 6610
Fax: (61-2) 9550 2841

Mindray Medical Nigeria

Mindray Medical India Pvt. Ltd. (Delhi Office)

16th floor, Building No. 9-B, DLF Cyber City, DLF Phase- 3,
Gurgaon- 122002, Haryana
Voice: (91-124) 4632-488
Fax: (91-124) 4632499

Mindray Medical India Pvt. Ltd. (Mumbai office)

B-404, City Point, Andheri Kurla Road, Andheri East, Mumbai,
400 059, India
Voice: (91-22) 4020-0000
Fax: (91-22) 4020-0011

Mindray Medikal Teknoloji Istanbul Ltd. Sti.

KOYALTI MEVKII 29 EKİM CAD.N.1 KUYUMCUKENT 2.NO.PLZ.3.KT.
D.1- 9 34520 BAĞÇELIEVLER/YENIBOSNA, İSTANBUL, TÜRKİYE
Voice: (90-212) 482 08 77
Fax: (90-212) 482 08 78

Mindray Medikal Teknoloji Istanbul Ltd. Sti. Branch Office Ankara

TURAN GÜNEŞ BULVARI ALEKSANDER DÜBÇEK CAD. NO:18/A
ÇANKAYA/ANKARA, TÜRKİYE
Voice: (90-312) 441 96 92
Fax: (90-312) 441 96 93

Mindray Medical Vietnam Co.,Ltd.


Unit 2, 4th Floor, "B" Building, Broadway Office Park,
102 Nguyen Luong Bang Blvd., Tan Phu Ward, Dist. 7, HCMC, Vietnam
Voice: (84) 8 5413-6970
Fax: (84) 8 5413-6971

Mindray Medical Thailand Limited

1768 Thai summit Tower Room 407, 4th floor, New petchaburee Road,
Kwaeng Bangkok, Khet Huai Khwang, Bangkok 10320, Thailand
Voice: (66) 2652-7288
Fax: (66) 2652-7266

Mindray Medical Korea

KNI Tower 9F 823-19 Yeoksam dong, Seoul-ganum, Korea
Tel: 008225688040
Fax: 008225688042

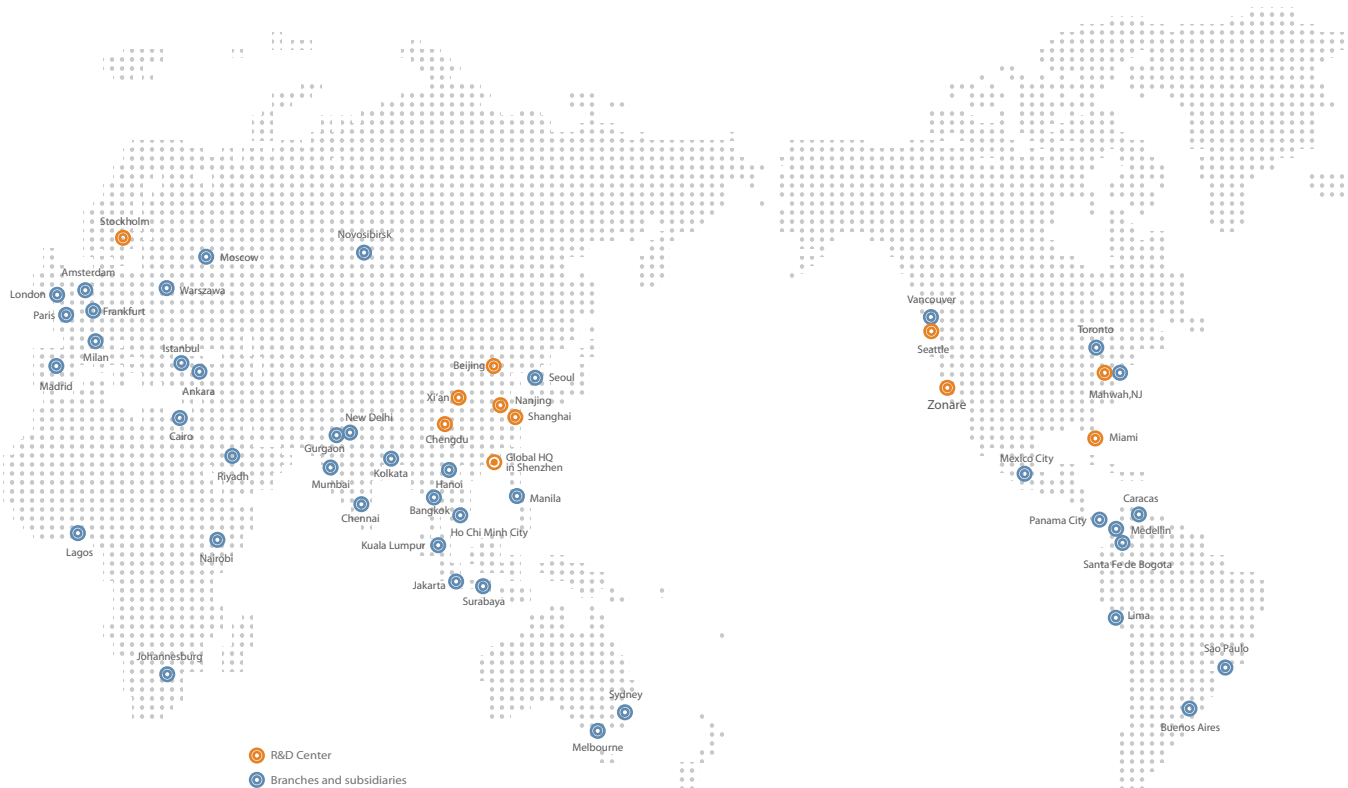

<p>Mindray Building, Keji 12th Road South, High-tech Industrial Park, Nanshan, Shenzhen 518057, P.R. China Tel: +86 755 8188 8998 Fax: +86 755 26582680 E-mail: intl-market@mindray.com www.mindray.com</p>	<p>mindray <i>healthcare with reach</i> are registered trademarks or trademarks owned by Shenzhen Mindray Bio-medical Electronics Co., Ltd. ©2014 Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. Specifications subject to changes without prior notice. Please refer to user manual for manufacturer information P/N:ENG- new WATO EX-35-210X285X8PX20160226</p> <p style="text-align: right;">mindray</p>
---	--